

**Franklin** ®

**PDR-3001**

# *Medical Book Card*



**POCKET PDR<sup>®</sup>**

**User's Guide**



## Quick Reference Guide

### To Use the Menus

1. Press **MENU** and highlight the Drugs menu or the Class menu.
2. Highlight an entry by typing its name or pressing **↑** or **↓**.
3. Press **ENTER**.
4. Select subitems, if needed, by highlighting them and pressing **ENTER**.
5. Use **↓**, **DN**, or **SPACE** to scroll.

### To Use Notes

1. When the text of the book is on screen, press **LIST**.
2. Press **ENTER** to select *Add note*.
3. Press **ENTER** or type your own note title and then press **ENTER**.
4. To find a note in the text, press **LIST**, highlight a note title, and then press **ENTER**.

### To Search for Words

1. Press **SEARCH**.
2. Type one or more words.
3. Press **ENTER**.
4. Hold **☆** and press **DN** to see the locations of matches.
5. Press **ENTER** to see a match or **MENU** to see the matches listed on the menus.
6. To expand a search, press **MORE**.

### Quick Keys\*

- (A)** Goes to *Adverse Reactions* topic.
- (C)** Goes to *Contraindications* topic.
- (D)** Goes to *Dosage and Administration* topic.
- (E)** Goes to *Equivalents* menu.
- (G)** Goes to *Generics* topic.
- (H)** Goes to *How Supplied* topic.
- (I)** Goes to *Indications and Usage* topic.
- (M)** Goes to *Manufacturer* topic.
- (S)** Goes to *How Supplied* topic.
- (T)** Goes to *Therapeutics Class Index*.
- (W)** Goes to *Warnings* topic.

\* You can use the Quick Keys only when the text of a drug monograph is on the screen.

## License Agreement

READ THIS LICENSE AGREEMENT BEFORE USING THE MEDICAL BOOK SYSTEM.

YOUR USE OF THE MEDICAL BOOK SYSTEM DEEMS THAT YOU ACCEPT THE TERMS OF THIS LICENSE. IF YOU DO NOT AGREE WITH THESE TERMS, YOU MAY RETURN THIS PACKAGE WITH PURCHASE RECEIPT TO THE DEALER FROM WHICH YOU PURCHASED THE MEDICAL BOOK SYSTEM AND YOUR PURCHASE PRICE WILL BE REFUNDED. The MEDICAL BOOK SYSTEM means the software product, hardware, and documentation found in this package and FRANKLIN means Franklin Electronic Publishers, Inc.

### LIMITED USE LICENSE

All rights in the MEDICAL BOOK SYSTEM remain the property of FRANKLIN. Through your purchase, FRANKLIN grants you a personal and nonexclusive license to use the MEDICAL BOOK SYSTEM. You may not make any copies of the MEDICAL BOOK SYSTEM or of the preprogrammed data stored therein, whether in electronic or print format. Such copying would be in violation of applicable copyright laws. Further, you may not modify, adapt, disassemble, decompile, translate, create derivative works of, or in any way reverse engineer the MEDICAL BOOK SYSTEM. You may not export or re-export, directly or indirectly, the MEDICAL BOOK SYSTEM without compliance with appropriate governmental regulations. The MEDICAL BOOK SYSTEM contains Franklin's confidential and propriety information which you agree to take adequate steps to protect from unauthorized disclosure or use. This license is effective until terminated. This license terminates immediately without notice from FRANKLIN if you fail to comply with any provision of this license.



## Contents

Key Guide .....	4
Installing Book Cards .....	5
Selecting a Book .....	6
Changing the Settings .....	6
Viewing a Demonstration .....	7
Finding Drug Information .....	7
Finding Therapeutic Classes .....	9
Searching for Words .....	10
Highlighting Search Words .....	11
Expanding Your Searches .....	12
Reading Special Markers .....	13
Using Notes .....	14
Looking Up Words in Other Books .....	15
Resetting the Medical Book System .....	15
Book Card Care .....	16
Specifications and Patents .....	16
Limited Warranty (U.S. only) .....	17
Index .....	18

### ► About Book Card Compatibility

**IMPORTANT** This book card can be used *only* with the Medical Book System™ platform and will not function in any other Franklin BOOKMAN® platform.

The Medical Book System platform can use Franklin BOOKMAN book cards as well as Medical Book System book cards. However, Franklin BOOKMAN platforms cannot use Medical Book System book cards.

### ► For More Information

To learn more about the Medical Book System or other products from Franklin Electronic Publishers, call 800-665-5450 or visit the Franklin Web site at [www.franklin.com](http://www.franklin.com).

## Publishers' Note

This Pocket PDR® book card contains excerpts of product information provided by drug manufacturers covering indications and usage, contraindications, warnings, adverse reactions, dosage and administration, and how supplied. Ophthalmology drugs are included in this 1998 Pocket PDR®. Precautionary material, information on clinical pharmacology, and certain other items are not included. Be sure to consult the printed edition of PDR®, PDR® Supplements, and PDR for Ophthalmology® in the event of questions.

Product information published in Pocket PDR® has been prepared by the manufacturers of the drugs. The function of the publisher is solely the compilation, organization, and distribution of this information. Drugs listed in Pocket PDR® are available only by prescription. Patients should consult a physician before using any of the medications.

In making manufacturers' information available, it should be understood that the publisher is not advocating the use of any product or dosage described. And, while every effort has been made to assure that this electronic version conforms with the drug manufacturers' originals, neither the pub-

lisher, Medical Economics Company, Inc., nor the producer of the unit, Franklin Electronic Publishers, Inc., warrant that the information is accurate, that the unit will perform without interruption or free of errors, or that the information contained therein satisfies government regulations requiring drug product information disclosure, nor shall they be liable for any claim arising out of or relating to use of Pocket PDR®.

For additional information on limitations of liability and warranty disclaimers, consult the License Agreement and Limited Warranty printed in this User's Guide. All use of this product is subject to the terms and conditions of the User License. Information in this book card is current through the 1998 editions of PDR® and PDR for Ophthalmology®. For later revisions, always check the latest PDR®, PDR® supplement, or PDR for Ophthalmology®.

Pocket PDR®.

Copyright © 1998 Medical Economics Company, Inc.

PDR® is a registered trademark used herein under license.

Copyright © 1998 Franklin Electronic Publishers, Inc.

All rights reserved, worldwide.



## Key Guide

### Color Keys

- MORE** (red) Expands a word search.
- SPEC** (green) Displays the Outline location of the current text.
- SEARCH** (yellow) Displays the word search screen.
- LIST** (blue) Lets you add, find, or remove a note.

### Function Keys

- BACK** Erases typed letters, or backs up to the previous screen.
- CAP** Shifts keys to type capitals or punctuation.
- CARD** Exits the currently selected book.
- CLEAR** Clears an entry or search and goes to the default state.
- ENTER** Enters a word, selects a menu item, or starts the highlight in text.
- HELP** Displays help messages.
- MENU** Displays the main menus.
- ON/OFF** Turns the platform on or off.
- SPACE** At entry screens, types a space. At menus and text, pages down.
- ? \*** At menus, displays the title of a highlighted item. With **CAP**, types an asterisk (\*) to stand for letters in a word. At text, displays the current Outline location of the text.

### Key Combinations\*

- ☆** + **BACK** Goes to the top level of a menu from a lower level.
- ☆** + **CARD** Transfers a highlighted word between installed book cards.
- ☆** + **ENTER** At a highlighted menu item, displays the related text. At text, highlights special markers (e.g., footnotes, cross-references, tables, etc.).
- ☆** + **Q-P** Types numbers.
- ☆** + **DN** or **UP** At text, displays the next or previous paragraph, table row, or search match.
- CAP** + **↑** or **↓** Goes to the top or bottom of a menu or list. At text, goes to the next or previous Outline title.

\* Hold the first key while pressing the second.

### ► Understanding the Keys

The functions of the keys may vary according to which book card is installed and selected in the platform. To learn how to use a particular book card, read its User's Guide.

## Key Guide

### Direction Keys

- UP** **DN** Pages up or down.
- ↕** Moves the cursor, text, or highlight.

### Quick Keys\*

- A** Goes to *Adverse Reactions* topic.
- C** Goes to *Contraindications* topic.
- D** Goes to *Dosage and Administration* topic.
- E** Goes to *Equivalents* menu.
- G** Goes to *Generics* topic.
- H** Goes to *How Supplied* topic.
- I** Goes to *Indications and Usage* topic.
- M** Goes to *Manufacturer* topic.
- S** Goes to *How Supplied* topic.
- T** Goes to *Therapeutics Class* index.
- W** Goes to *Warnings* topic.

\* You can use the Quick Keys only when the text of a drug monograph is on the screen.


### ► Understanding the Quick Keys

When a drug monograph is on the screen, you can use the quick keys to quickly view its different parts, its equivalents, or the Therapeutics Class index.

## Installing Book Cards

**CAUTION** Never install or remove a book card when the platform is turned on. If you do, information entered in any installed book cards may be erased.

1. Turn the platform off.
2. Turn the platform over.
3. Align the book card tabs with the notches in a card slot.

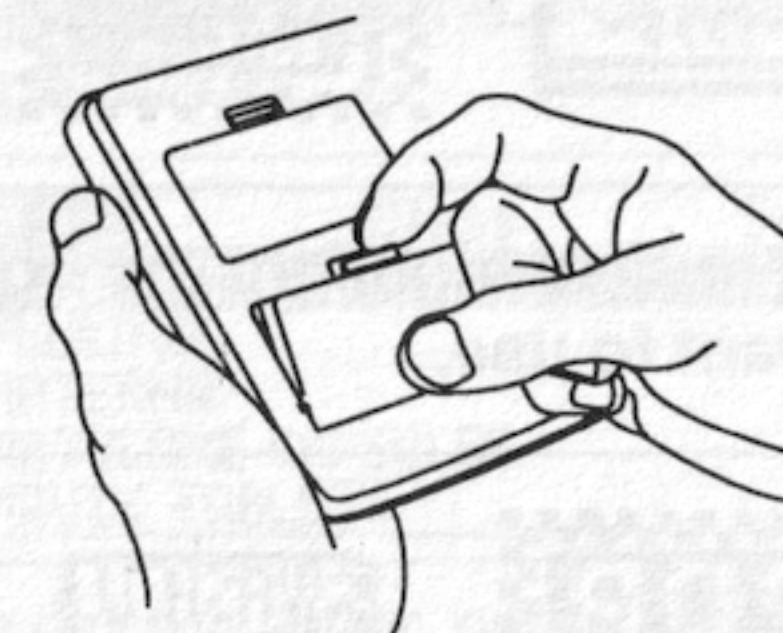

4. Press the book card until it snaps into the slot.

### ► Removing Book Cards

**CAUTION** When you remove a book card from a platform, information entered in that book card may be erased.



## Selecting a Book

If you have installed two book cards in the platform, you can select which book you want to read.

1. Turn the platform on.

2. Press **CARD**.



3. Use **←** or **→** to highlight the book you want to use.



4. Press **ENTER** to select it.

### ► Resuming Where You Left Off

You can turn off the platform at any screen. When you turn it on again, the last screen that you viewed appears.

## Changing the Settings


Using the Setup menu, you can adjust the type size, default topic, default state, shutoff time, and screen contrast of this book.

The type size sets how large the characters appear on screen. The shutoff time sets how long your Medical Book System stays on if you forget to turn it off.

The default topic sets which topic will appear first when you hold **☆** and press **ENTER** to view a drug monograph. The default state sets the screen that appears when you press **CLEAR**.

1. Press **MENU**.

2. Highlight *Set Type Size*, *Set Default Topic*, *Set Default State*, *Set Shutoff*, or *Set Contrast* on the Setup menu.

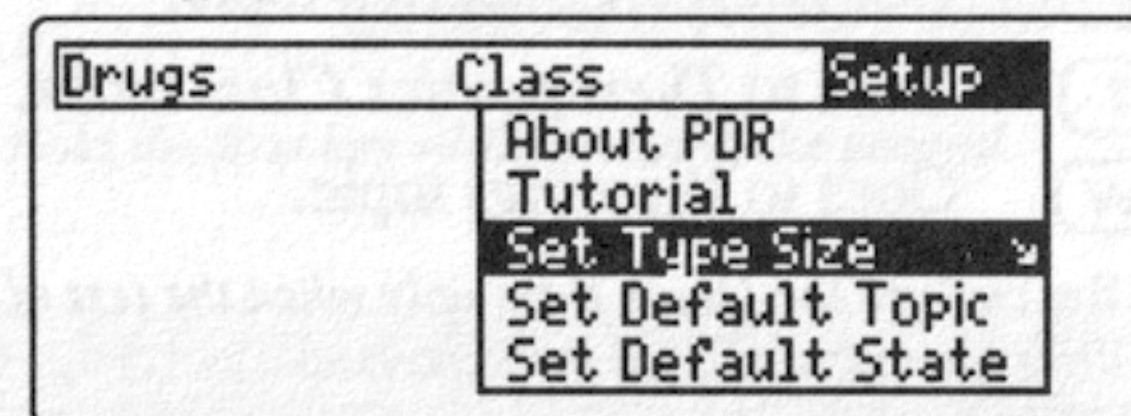

3. Press **ENTER**.

4. Use **↓** or **↑** to change the setting.


Or press **BACK** to leave the setting unchanged.

5. Press **ENTER** to select it.

## Viewing a Demonstration

Before you start using this book, you may want to see a brief demonstration.

Press **MENU** and then use the arrow keys to highlight *View Demo* on the Setup menu.

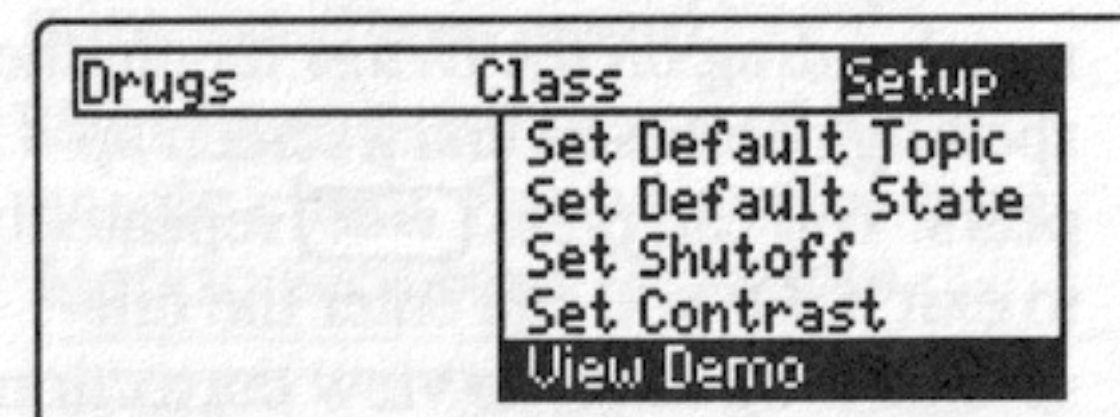

Press **ENTER** to select it. To stop the demonstration, press **CLEAR**.

### ► Help is Always at Hand

You can view a help message at virtually any screen by pressing **HELP**. To exit help, press **BACK**.

To read a tutorial about this book, select *Tutorial* from the Setup menu.

For more information about the Pocket PDR®, select *About PDR* from the Setup menu.

### ► Follow the Arrows

The flashing arrows at the right of the screen show which arrow keys you can press to move through menus or view more text.


## Finding Drug Information

The simplest way to find drug information is to type a brand or generic drug name at the Drugs menu. Brand names are capitalized on the Drugs menu, but you do not need to type capitals.

1. Press **CLEAR**.

2. If needed, press **MENU** and highlight the Drugs menu.

3. Type a drug name. For example, type *Marax*.

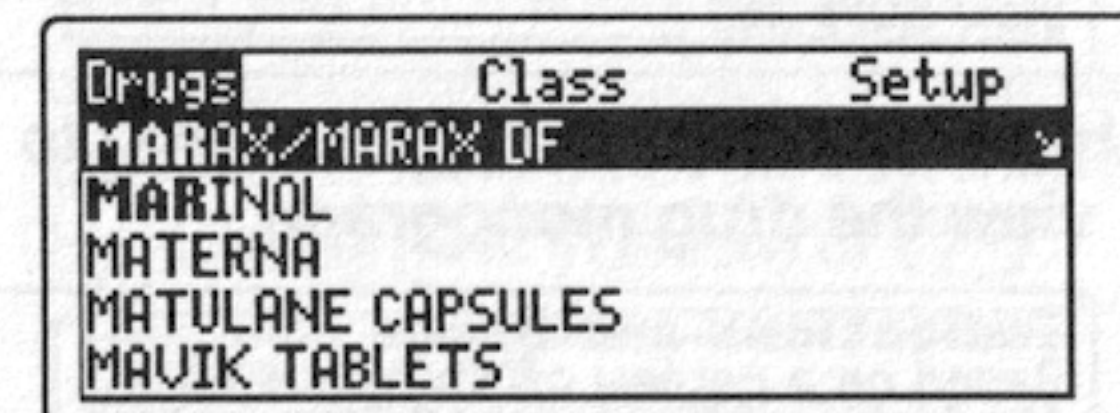

To delete a letter, press **BACK**.

You can also highlight a menu item by using **↓** or **↑**.

The slanted arrow indicates that the menu item has a submenu.

4. When the drug is highlighted, press **ENTER** to select it.

If you entered a generic name, select one of the brand name drugs that contain the selected generic name.



## Finding Drug Information

If you entered a brand name, the monograph's topics appear.

Note: You can go directly from the Drugs menu to the default topic of the monograph by holding **[☆]** and pressing **[ENTER]**.

MARAX/MARAX DF  
Indications and Usage  
Contraindications  
Adverse Reactions  
Dosage and Administration  
How Supplied

5. Highlight a topic and press **[ENTER]** to view the drug monograph.

**Indications and Usage**  
Based on a review of this drug by the National Academy of Sciences-National Research Council and/or other information, FDA has classified the indications as follows:  
"Possibly" Effective: For controlling bronchospastic disorders.

6. To go quickly to another topic, use the Quick Keys.

To learn more about the Quick Keys, read the "Key Guide."

7. Use **[↓]**, **[DN]**, or **[SPACE]** to scroll.  
8. To go to the next or previous topic, hold **[CAP]** and use **[↓]** or **[↑]**.

## Finding Drug Information

9. Press **[CLEAR]** when finished.

### ► If You Misspell a Drug Name

When you type letters that do not match a drug on the *Drugs* menu, the spelling correction entry screen appears. You can press **[BACK]** repeatedly to exit it. Or you can enter the misspelled drug name to view corrections.

Corrections for digabind  
Words starting with digabind  
DIGABIND  
dioxybenzone

You can select *Words starting with...* to view completions, if any, or select a correction to view its monograph.

### ► Finding Drug Equivalents and Therapeutic Classes

When reading a drug monograph, you can quickly find its equivalents or a list of drugs in its therapeutic class. To find a drug's equivalents, if any, press **[E]** and then select a monograph. To find a list of drugs in its therapeutic class, press **[T]**.

## Finding Therapeutic Classes

Every drug in the Pocket PDR® is classified by its therapeutic use.

1. Press **[CLEAR]**.
2. If needed, press **[MENU]** and highlight the Class menu.
3. Start typing the name of a therapeutic class. For example, type *Motion Sickness Products*.

Drugs	Class	Setup
MOTION SICKNESS PRODUCTS		
MULTIPLE SCLEROSIS MANAGEMENT		
MUSCLE RELAXANTS		
NARCOTIC DETOXIFICATION		
NASAL PREPARATIONS		

To delete a letter, use **[BACK]**.

You can also use **[↓]** or **[↑]** to highlight a class.

4. When the class is highlighted, press **[ENTER]**.
5. Highlight a drug and press **[ENTER]**. For example, select *Phenergan Syrup*.
6. Highlight a topic and press **[ENTER]**.

**Contraindications**  
Promethazine is contraindicated in individuals known to be hypersensitive or to have had an idiosyncratic reaction to promethazine or to other phenothiazines.  
Antihistamines are contraindicated for use in the treatment of lower

## Finding Therapeutic Classes

7. To go quickly to another topic, use the Quick Keys.

To learn more about the Quick Keys, read the "Key Guide."

8. Use **[↓]**, **[DN]**, or **[SPACE]** to scroll.
9. To go to the next or previous topic, hold **[CAP]** and use **[↓]** or **[↑]**.
10. Press **[CLEAR]** when finished.

### ► Where in the Book Am I?

When you are reading the text of a drug monograph or table, you can see the Outline location of the current text by pressing **[SPEC]** (green) or **[?]**.

Outline location:  
MARAX/MARAX DF  
Indications and Usage

To scroll up to a higher level of the Outline, use **[↑]**.

To display the full title of an item, press **[?]**. This is especially useful when the title ends in ellipses.

Press **[BACK]** to go back to the text that you were reading. Or press **[ENTER]** to go to the text of the highlighted item.



## Searching for Words

Word searches find all the occurrences, or matches, of a word or word phrase in this book.

You can type up to 99 characters in a single search. But you cannot search for short, common words such as *the* and *is*.

1. Press **SEARCH** (yellow).
2. Type your search words. For example, type *flu*.

Type words, press **ENTER**.

flu

Press **MENU** for other options.

To delete a letter, use **BACK**.

You do not need to type capitals or punctuation to find matches.

3. Press **ENTER** to search.

Outline location:

ALLEGRA  
Adverse Reactions  
**TABLE Adverse Experiences...**  
Viral Infection (cold, flu)

This is the Outline location of the first match of your search words.

4. Press **ENTER** to see the match.

## Searching for Words

- **Viral Infection (cold, flu)**
  - Fexofenadine 60 Mg Twice Daily (n=679): 2.5%
  - Placebo Twice Daily (n=671): 1.5%
- **Nausea**
  - Fexofenadine 60 Mg Twice Daily (n=679): 1.6%

Your search words are boxed.

5. To go quickly to another topic, use the Quick Keys.

To learn more about the Quick Keys, read the "Key Guide."

6. To see the Outline location of the next or previous matches, if any, hold **☆** and press **DN** or **UP**.

Outline location:

ALUPENT  
Adverse Reactions  
**TABLE ALUPENT® Tablets: I...**  
OTHER  
Flu symptoms

7. Press **MENU** to see the matches listed on the menus.

After a word search, only matching entries are listed on the menus.

8. Use the arrow keys to highlight a menu item and then press **ENTER**. If needed, select a subitem.
9. Press **SEARCH** (yellow) or **CLEAR** to clear your search.

## Searching for Words

### ► If You Misspell a Search Word

If you enter a misspelled search word, a list of corrections appears.

Corrections for **psoriasis**

psoriasis  
psoriasi  
serositis  
scissors  
psoriatic

Use **↓** to highlight a correction and then press **ENTER**.

### ► Searching for Parts of Words

If you want to find prefixes, suffixes, or other parts of words, type an asterisk (\*) in place of the missing letters. For example, type *bacter\** and then press **ENTER**. To type an asterisk, hold **CAP** and press **?\***.

Matches for **bacter\***

Bacteroides  
bacteremia  
bacteremias  
bacteria  
bacterial

Use **↓** to highlight a match and then press **ENTER**.

## Highlighting Search Words

You can also search for a word by highlighting it in text.

1. At the text, press **ENTER**.

**Contraindications**  
ANTI-KALIURETIC THERAPY AND POTASSIUM SUPPLEMENTATION  
Dyazide should not be given to patients receiving other potassium-sparing agents such as spironolactone, amiloride or other formulations containing triamterene. Con-

To turn off the highlight, press **BACK**.

2. Use the arrow keys to highlight a word.

**Contraindications**  
ANTI-KALIURETIC THERAPY AND POTASSIUM SUPPLEMENTATION  
Dyazide should not be given to patients receiving other potassium-sparing agents such as spironolactone, amiloride or other formulations containing **triamterene**. Con-

3. Press **ENTER** to search for it.
4. Press **ENTER** again to see the first match.
5. Hold **☆** and press **DN** repeatedly to see the Outline locations of the next matches, if any.  
To see the location of previous matches, if any, hold **☆** and use **UP**.
6. Press **MENU** to see the menu items with matches.
7. Press **CLEAR** when finished.



## Expanding Your Searches

You can expand a word search up to three times to find inflected forms of words, synonyms, and similar words.

1. Press **SEARCH** (yellow).
2. Type a word and then press **ENTER**.  
For example, enter *diabetes*.

Type words, press **ENTER**.

diabetes 4 Matches

Press **MENU** for other options.

3. Press **MORE** (red) to expand the search.

Outline location:

PROGRAF

Warnings

298 Matches

4. Hold **☆** and use **DN** to see the Outline locations of the next matches, if any.  
To see the location of previous matches, if any, hold **☆** and use **UP**.
5. Press **MENU** to see the matches listed on the menu.
6. Press **CLEAR** when finished.

## Expanding Your Searches

### ► Understanding Word Searches and Expansions

When you search for a word, either by entering it at the word search screen or by highlighting it in text, headings and subheadings (in bold type) are searched first for exact matches. If no exact matches are found in the headings, then the text is searched for exact matches.

If no exact matches are found, then the headings are searched for inexact matches, such as inflections and synonyms. If no inexact matches are found in the headings, then the text is searched for them.

A word search stops as soon as matches are found at a level described above. Pressing **MORE** (red) continues the search to the next level and finds additional matches, if any. Note: The number of matches found at a given search level may be less than at the previous level.

## Reading Special Markers

Throughout this book, you will find cross-references (indicated by **XREF**'s), footnotes (indicated by **FOOT**'s), bibliographic references (indicated by **BIB**'s), and tables (indicated by **TABLE**'s). Cross-references take you directly to related material in the book. Tables, footnotes, and bibliographic references provide additional information.

1. When you see a **XREF**, **FOOT**, **BIB**, or **TABLE** in the text, hold **☆** and press **ENTER**.

has been inadequate (see National Cholesterol Education Program (NCEP) Guidelines, summarized in Table 3).  
See table 3.  
At the time of hospitalization for an acute coronary event, consideration can be given to initiating drug

Notice the highlight. To undo the highlight, press **BACK**.

2. Press **ENTER**.

**TABLE Table 3: NCEP Guidelines for Lipid Management**

- **Definite Atherosclerotic Disease** **FOOT**: No
- Two or More Other Risk Factors **FOOT**: No
- LDL-Cholesterol:
  - Initiation Level:  $\geq 190$  mg/

3. Press **BACK** repeatedly to go back.

## Reading Special Markers

### ► Reading Multiple Special Markers

Sometimes more than one special marker appears on screen at one time.

To highlight the first special marker, hold **☆** and press **ENTER**. To highlight subsequent cross-references, footnotes, bibliographic references, and tables, continue holding **☆** and press **ENTER** again.

When you have highlighted the special marker that you want, press **ENTER** to read it. Then press **BACK** to return to the highlighted item.

### ► Using and Understanding Tables

In this book card, tables are displayed as bulleted text. Row titles are in bold type and column titles and information are indented beneath them.

When selected, most cross-reference markers (**XREF**) take you to a table. Once the table is in view, you can usually return to its cross-reference by highlighting and selecting the table marker (**TABLE**).



## Using Notes

### To Add a Note

You can add up to 10 notes to the text of this book to quickly find the sections that you most often use.

1. When text is on screen, press **LIST** (blue).
2. Press **ENTER** to select *Add note*.

Type your note then press **ENTER**. Or press **HELP**.

YABOKINASE

Press **STAR** **BACK** to cancel.

3. Do one of the following:
  - To accept the note title, press **ENTER**.
  - To change the note title, type your changes and then press **ENTER**.

Note titles can contain up to 24 characters.

- To cancel, hold **☆** and press **BACK**.

**NOTE** **YABOKINASE**  
**Indications and Usage**  
**PULMONARY EMBOLISM**  
YABOKINASE (urokinase for injection) is indicated in adults:  
-- For the lysis of acute massive pulmonary emboli, defined as obstruction of blood flow to a lobe

The **NOTE** marker is added to the text when you add a note.

## Using Notes

### To Find a Note

Once you have added notes to this book, you can easily find them.

1. Press **LIST** (blue).
2. Use **↓** to highlight a note title.  
Your notes are listed by recency.
3. Press **ENTER** to go to that note.

### To Remove a Note

1. Find a **NOTE** in the text of this book.  
To learn how, read "To Find a Note" above.
2. Press **LIST** (blue).
3. To see the full title of the note to be removed, press **?\***.
4. Press **ENTER** to select *Remove note...*  
Or press **BACK** to exit without removing the note.

### ► Removing All Your Notes

You can remove all your notes at one time by resetting the Medical Book System. To learn how, read "Resetting the Medical Book System."

**WARNING** Resetting the Medical Book System may erase other information that you have entered in other installed book cards.

## Looking Up Words in Other Books

The Medical Book System platform enables you to look up a word from one installed book card in the other installed book card.

This book card can send words to and receive words from other book cards. Some book cards are not able to transfer words. To learn if a book card can send or receive words, read its User's Guide.

1. Install two book cards in the platform.
2. Select a book to read.
3. Highlight a word in that book.

You can highlight words in menus, lists, entries, and other text. To learn how to highlight a word in this book, read "Highlighting Search Words."

4. Hold **☆** and press **CARD**.
5. Highlight the icon of the other book card and then press **ENTER**.  
The word that you highlighted appears in the other book.
6. Press **ENTER** again, if needed, to search for the word in the other book.

## Resetting the Medical Book System

If, due to electrostatic discharge or other cause, the keyboard fails to respond or the screen performs erratically, try the steps below to reset the platform. Perform only as many steps as required to restore normal operation.

**CAUTION** Resetting the platform when book cards are installed may erase settings and information entered in those book cards. Remove all book cards before resetting.

1. Hold **CLEAR** and press **ON/OFF** twice.  
If nothing happens, try Step 2.
2. Use the end of an opened paper clip to gently press the reset button on the back of the platform.

The reset button is recessed in a pin-sized hole located near the book card slots.

**CAUTION** Pressing the reset button with more than light pressure may permanently disable it.

If nothing happens, try Step 3.

3. Remove and reinstall the batteries.  
To learn how, read the Medical Book System platform's User's Guide.




## Book Card Care

- **Do not touch the metal contacts on the book cards.**

**CAUTION** Touching the electrical contacts with statically charged objects, including your fingers, could erase information entered in a book card.



- **Do not put excessive pressure on the book cards.**

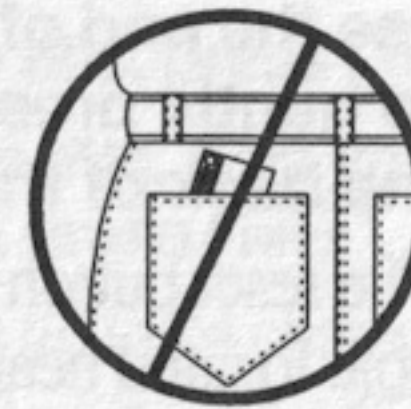

- **Do not expose the book cards to heat, cold, or liquids.**



## Specifications and Patents

### Model PDR-3001

- size: 13.6 x 9.2 x 2 cm
- weight: 0.4 oz

© 1998 Franklin Electronic Publishers, Inc. All rights reserved.

© 1998 Medical Economics Company, Inc. All rights reserved.

Physician's Desk Reference is a registered trademark used herein under license.

Medical Book System and BOOKMAN are trademarks of Franklin Electronic Publishers, Inc.

**FCC Notice:** Complies with the limits for a Class B computing device pursuant to Subpart B of part 15 of FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference; and (2) This device must accept any interference received, including interference that may cause undesired operation.

U.S. Patents 4,490,811; 4,830,618;  
5,113,340; 5,218,536; 5,321,609;  
5,396,606; 5,627,726;

Euro. Patent 0 136 379.

German Pats. M9409743.7, M9409744.5.  
PATENTS PENDING.

ISBN 1-56712-418-6

## Limited Warranty (U.S. only)

Franklin Electronic Publishers, Inc. ("Franklin") warrants to the end user that this product will be free from defects in material and workmanship for a period of ONE YEAR from the date of original retail purchase, as evidenced by sales receipt. On discovery of a defect, the end user must return this product (transportation charges prepaid) either to the dealer from whom it was purchased or directly to Franklin at the address given below. Each product returned must include the user's name, address, and telephone number, as well as a brief description of the nature of the defect and a copy of the sales receipt as proof of the date of the original retail purchase. Franklin will, at its option, repair or replace any product at no further charge to the end user on determination by Franklin, in its sole discretion, that the product was defective and that such defect arose within the duration of this limited warranty.

Enclose \$3.00 check or money order for shipping and handling to Franklin Electronic Publishers, Inc.

This warranty does not apply if, in the sole discretion of Franklin, the product has been tampered with, damaged by accident, abuse, misuse, or misapplication, or as a result of service or modification by any party, including any dealer, other than Franklin. FRANKLIN SHALL NOT BE RESPONSIBLE IN ANY WAY FOR ANY SERVICE OR MODIFICATION TO THIS PRODUCT BY ANY PARTY, INCLUDING ANY DEALER, OTHER THAN FRANKLIN.

This warranty applies only to products manufactured by or for Franklin. Batteries, corrosion of battery contacts and any damage caused by batteries are not covered by this warranty.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, SHALL BE STRICTLY LIMITED IN DURATION TO THAT OF THE EXPRESS WARRANTY SET FORTH ABOVE, THAT IS, ONE YEAR FROM THE DATE OF ORIGINAL RETAIL PURCHASE.

THE WARRANTY AND REMEDY SET FORTH ABOVE ARE THE EXCLUSIVE REMEDY OF THE RETAIL

BUYER AND END USER IN CONNECTION WITH THE MANUFACTURE, SALE, OR USE OF THIS PRODUCT AND ARE IN LIEU OF ANY AND ALL OTHER WARRANTIES OR REMEDIES, WRITTEN OR ORAL, EXPRESS OR IMPLIED. NO FRANKLIN DEALER, AGENT, OR EMPLOYEE IS AUTHORIZED TO MAKE ANY ADDITIONAL WARRANTY IN THIS REGARD OR TO MAKE ANY MODIFICATION OR EXTENSION OF THIS EXPRESS WARRANTY.

FRANKLIN SHALL NOT BE RESPONSIBLE FOR SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM THE BREACH OF ANY EXPRESS OR IMPLIED WARRANTY OR IN CONNECTION WITH THE MANUFACTURE, SALE, OR USE OF THIS PRODUCT UNDER ANY LEGAL THEORY, WHICH DAMAGES SHALL INCLUDE, BUT NOT BE LIMITED TO, LOST PROFITS, DAMAGES TO PROPERTY OR DAMAGES FOR PERSONAL INJURY (BUT ONLY TO THE EXTENT PERMITTED BY LAW). UNLESS OTHERWISE CONTRARY TO APPLICABLE LAW, FRANKLIN'S LIABILITY SHALL IN NO CASE EXCEED THE PRICE PAID FOR THE PRODUCT CLAIMED TO BE DEFECTIVE.

This warranty shall not be applicable to the extent that the enforcement of any provision may be prohibited by applicable law. This warranty gives you specific rights, and you may also have other rights which vary from state to state. Some states do not allow the exclusion or limitation of incidental or consequential damages or limitation on how long an implied warranty may last so the above limitation or exclusion may not apply to you. Enclose \$3.00 check or money order for shipping and handling to:

FRANKLIN ELECTRONIC PUBLISHERS, INC.  
ATTENTION: SERVICE DEPARTMENT  
ONE FRANKLIN PLAZA  
BURLINGTON, NEW JERSEY 08016-4907  
609-386-2500



This product, excluding batteries, is guaranteed by Franklin for a period of one year from the date of purchase. It will be repaired or replaced (at Franklin's option) free of charge for any defect due to faulty workmanship or materials.

Products purchased outside the United States that are returned under warranty should be returned to the original vendor with proof of purchase and description of fault. Charges will be made for all repairs unless valid proof of purchase is provided.

This warranty explicitly excludes defects due to misuse, accidental damage, or wear and tear. This guarantee does not affect the consumer's statutory rights.

This unit may change operating modes due to Electrostatic Discharge. Normal operation of this unit can be re-established by pressing the reset key, **ON/OFF**, or by removing/replacing batteries.

**?\* key** 4, 9, 11, 14  
**Arrow keys** 5, 7  
**Auto-resume feature** 6  
**Automatic shutoff** 6  
**BIB marker** 13  
**Book cards**  
     installing 5  
     platform compatibility 2  
     removing 5  
     selecting 6  
**Changing**  
     screen contrast 6  
     shutoff time 6  
     type size 6  
**Color keys**  
     key guide 4  
     understanding 4  
**Default state**  
     changing 6  
     defined 6  
**Demonstration** 7  
**Direction keys** 5  
**Equivalents, finding** 8  
**FOOT marker** 13  
**Function keys** 4  
**Green (SPEC) key** 9  
**Help messages** 7

**Highlighting**  
     bibliographic references 13  
     footnotes 13  
     words 11  
**Key combinations** 4  
**Misspellings, correcting**  
     drug names 8  
     search words 11  
**MORE (red) key** 12  
**NOTE marker** 14  
**Notes**  
     deleting 14  
     deleting all at once 14  
     finding 14  
**Pocket PDR®**  
     description of 3  
**Quick Keys** 5, 8, 9  
**Red (MORE) key** 12  
**Reset button** 15  
**Screen contrast** 6  
**SEARCH (yellow) key** 10, 12  
**Shutoff time** 6  
**SPEC (green) key** 9  
**Special markers, reading** 13  
**TABLE marker** 13  
**Tables**  
     understanding 13  
     using 13  
**Therapeutic classes, finding** 8, 9  
**Tutorial** 7  
**Type size** 6  
**Viewing**  
     bibliographic references 13  
     cross-references 13  
     demonstration 7  
     help messages 7  
**Word searches, understanding** 12  
**Words, highlighting** 11  
**XREF marker** 13  
**Yellow (SEARCH) key** 10, 12



**Franklin®**



*Medical Book System™*